


PARMAX
D E N T A L P R O D U K T E R

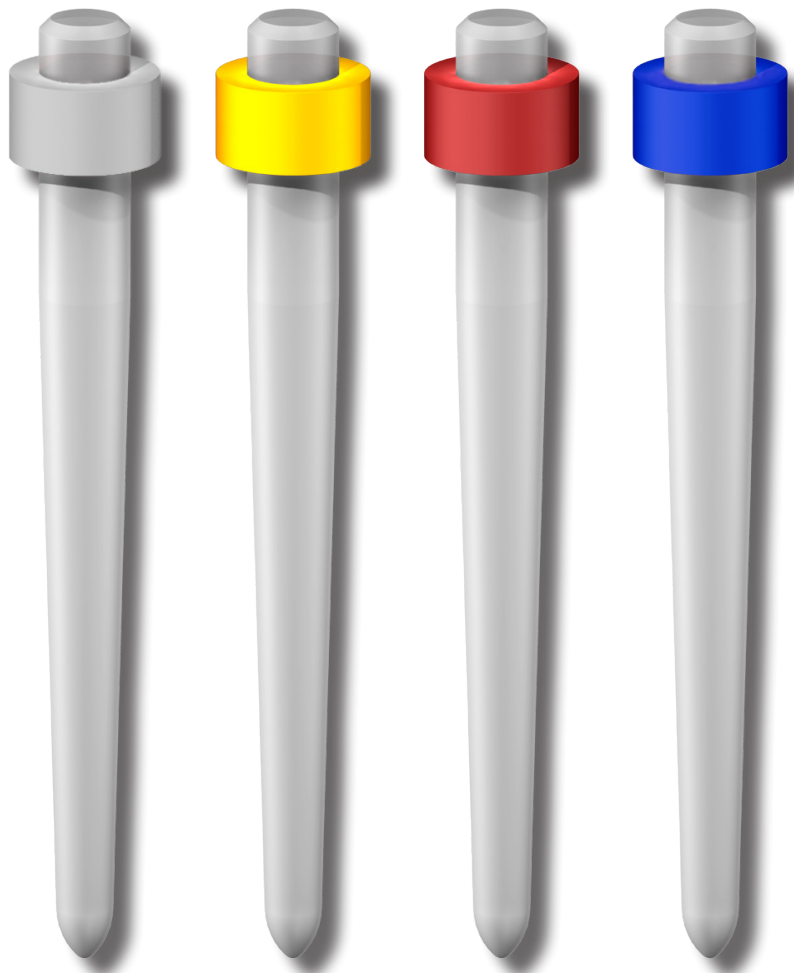

Product Catalog 2025



TABLE OF CONTENTS

TERMS & CONDITIONS	3
-------------------------------	----------

PARMAX POSTS	4
---------------------	----------

Metal Posts

Titanium	5
Gold Plated	6
Reamers & Keys	7

Fiber Posts

Paralight Posts	9
Paralight Reamers	9

PARMAX DIAMONDS	10
------------------------	-----------

Parmax Classic Multi-Use Diamonds	10
Parmax Polishing Diamonds	11
Parmax Single-Use Diamonds	11



TERMS & CONDITIONS

PRICES

Prices are subject to change without notice. Possession of a price list does not constitute an offer to sell. All prices are in European currency, Euro (€). For sales within the EU, Swedish VAT (25%) will be added unless a valid VAT registration number is presented at ordering.

ORDERING

When ordering products, please include complete catalog numbers.

MINIMUM ORDERS

Orders totaling less than € 50 will be subject to a special processing charge of € 20.

TERMS

Payment will be made against documents unless otherwise agreed upon. Past due balances are subject to interest according to EU regulations. Dentatus also reserves the right to renegotiate the terms and conditions in the event of repeated late payments. Delivered products remain under the property of Dentatus until complete payment has been made by the customer.

FREIGHT/SHIPPING

GENERAL

All shipments are Ex Works (EXW). Assistance with arranging freight and insurance on request can be provided.

DOCUMENTS

Costs for documents required to be authenticated and/or legalized by national authorizations will be invoiced separately for each shipment.

CLAIMS

Short claims must be made within 10 days from receipt of goods. Damaged goods must be reported to the carrier within 24 hours. Damage claim must be filled by consignee.

Parmax Metal Posts

The Root Shaped Threaded Posts

High Quality Metal Posts with Anatomical Geometry for Passive Placement


The Parmax Metal Posts have an anatomical geometry for passive placement and cementation. The shape is recognized by many as the standard in dental practice. Parmax Metal Posts are available in biocompatible titanium and economical gold plated brass. Its head retention can be customized - the square cross-cut head can be expanded in any plane to fit within the anticipated aesthetic core form.

PARMAX CLASSIC POST SYSTEM

Parmax Metal Posts are based on the shape of the classic root shaped threaded posts. This is an economical alternative for a permanent solution. Available in assortments and refills.

SIMPLE TECHNIQUE

1. The canal is explored with a small instrument, leaving a 4-5 mm apical sealed end.
2. A selected Parmax reamer size is used to widen the canal to intended depth.
3. The Metal Post matching the reamer size is tested and may be angled to the desired position for the core and crown.
4. The post is luted in place with a resin material for a built-up core for cemented crowns.

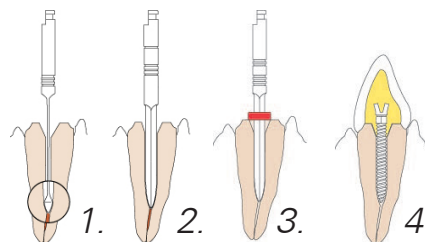

			1	2	3	4	5	6
Ø mm			1,05	1,20	1,35	1,50	1,65	1,80
pcs.								
S	7.8 mm	12	PT-S1	PT-S2	PT-S3	PT-S4	PT-S5	PT-S6
M	9.3 mm	12	PT-M1	PT-M2	PT-M3	PT-M4	PT-M5	PT-M6
L	11.8 mm	12	PT-L1	PT-L2	PT-L3	PT-L4	PT-L5	PT-L6
EL	14.2 mm	12			PT-EL3	PT-EL4	PT-EL5	PT-EL6
SL	17.0 mm	12			PT-SL3	PT-SL4	PT-SL5	PT-SL6

PARMAX METAL POSTS - TITANIUM

Art no:

PT-60

Product name:
Starter kit, no reamers

60 anchors (20 sizes, 3 of each)

1 hollow key, 1 cross key, 20 stoppers



Parmax Titanium Starter Kit

Refills, 12 pcs/pkg

Short, 7,8 mm

PT-S1	Diameter, 1.05 mm
PT-S2	Diameter, 1.20 mm
PT-S3	Diameter, 1.35 mm
PT-S4	Diameter, 1.50 mm
PT-S5	Diameter, 1.65 mm
PT-S6	Diameter, 1.80 mm

Medium, 9,3 mm

PT-M1	Diameter, 1.05 mm
PT-M2	Diameter, 1.20 mm
PT-M3	Diameter, 1.35 mm
PT-M4	Diameter, 1.50 mm
PT-M5	Diameter, 1.65 mm
PT-M6	Diameter, 1.80 mm

Long, 11,8 mm

PT-L1	Diameter, 1.05 mm
PT-L2	Diameter, 1.20 mm
PT-L3	Diameter, 1.35 mm
PT-L4	Diameter, 1.50 mm
PT-L5	Diameter, 1.65 mm
PT-L6	Diameter, 1.80 mm

Extra long, 14,2 mm

PT-EL3	Diameter, 1.35 mm
PT-EL4	Diameter, 1.50 mm
PT-EL5	Diameter, 1.65 mm
PT-EL6	Diameter, 1.80 mm



Parmax Titanium

PARMAX METAL POSTS - GOLD PLATED

Art no:
Product name:

PG-60

Starter kit, no reamers

60 anchors (20 sizes, 3 of each)

1 hollow key, 1 cross key, 20 stoppers

PG-240

Starter kit, no reamers

240 anchors (20 sizes, 12 of each)

1 hollow key, 1 cross key, 20 stoppers



Parmax Gold Plated Starter Kit

Refills, 12 pcs/pkg

Short, 7,8 mm

PG-S1

Diameter, 1.05 mm

PG-S2

Diameter, 1.20 mm

PG-S3

Diameter, 1.35 mm

PG-S4

Diameter, 1.50 mm

PG-S5

Diameter, 1.65 mm

PG-S6

Diameter, 1.80 mm

Medium, 9,3 mm

PG-M1

Diameter, 1.05 mm

PG-M2

Diameter, 1.20 mm

PG-M3

Diameter, 1.35 mm

PG-M4

Diameter, 1.50 mm

PG-M5

Diameter, 1.65 mm

PG-M6

Diameter, 1.80 mm

Long, 11,8 mm

PG-L1

Diameter, 1.05 mm

PG-L2

Diameter, 1.20 mm

PG-L3

Diameter, 1.35 mm

PG-L4

Diameter, 1.50 mm

PG-L5

Diameter, 1.65 mm

PG-L6

Diameter, 1.80 mm

Extra long, 14,2 mm

PG-EL3

Diameter, 1.35 mm

PG-EL4

Diameter, 1.50 mm

PG-EL5

Diameter, 1.65 mm

PG-EL6

Diameter, 1.80 mm



Parmax Gold Plated

PARMAX CLASSIC REAMERS

Art no:
Product name:

PRA-1/6

Assortment, short reamers

6 pcs, one of each size

Refills, 3 pcs/pkg

Short, 28 mm

PRA-1/3FP

Diameter, 1.05 mm

PRA-2/3FP

Diameter, 1.20 mm

PRA-3/3FP

Diameter, 1.35 mm

PRA-4/3FP

Diameter, 1.50 mm

PRA-5/3FP

Diameter, 1.65 mm

PRA-6/3FP

Diameter, 1.80 mm


Reamers, short

PRB-1/6

Assortment, long reamers

6 pcs, one of each size

Refills, 3 pcs/pkg

Long, 33 mm

PRB-1/3FP

Diameter, 1.05 mm

PRB-2/3FP

Diameter, 1.20 mm

PRB-3/3FP

Diameter, 1.35 mm

PRB-4/3FP

Diameter, 1.50 mm

PRB-5/3FP

Diameter, 1.65 mm

PRB-6/3FP

Diameter, 1.80 mm


Reamers, long

ACCESSORIES - KEYS

Art no:
Product name:

P-2K

Cross Key, 2 pcs


P-2K

P-2H

Hollow Key, 2 pcs


P-2H

Parmax Paralight Posts

Fiberglass Resin Post and Core System



For Strong Esthetic Restorations

Parmax Paralight posts offer the One-Step procedure for simultaneously curing of adhesives and composites within the confines of canals, creating a mono-bloc strength anchor-core support.

The resin-reinforced Fiberglass Paralight posts are the optimal choice for anchorage and a cohesive foundation for perfect esthetics. The anchors can be economically placed in the smallest and largest canals with off-the-shelf brand name adhesives and resins. Paralight transmits natural hues and tooth colors through the teeth and the surrounding gingival tissue.

ESTHETICS

The Parallax Paralight Posts transmit light into the confines of canals, polymerizing adhesives and composites attaching firmly the anchor with core foundation. The Paralight Posts, with modern adhesive materials, overcome the disadvantages of metal post products especially in the highly esthetic visible zones. It creates the ideal color background for crowns and the surrounding gingival tissue that is accomplished with conventional techniques and procedures.

ADHESION

The Paralight Posts are composed of longitudinal fiberglass rods in a strong resin matrix with a modulus of elasticity well within the range of dentin. It provides the strength needed for supporting restorations. The anchor-core construction creates the ideal mono-bloc strength foundation for aesthetic restorations.

TECHNICAL INFORMATION

Paralight provide extremely high adhesion to composite materials. The post and core is retentive and antirotational, conserving sound radicular and coronal structure. Available in 3 sizes, S (Ø 1,45 mm), M (Ø 1,6 mm), L (Ø 1,75 mm).

PARMAX PARALIGHT FIBERGLASS POSTS

Art no:	Product name:
PL-XS15	Paralight 15 posts, x-small, with Depth Marker
PL-S15	Paralight 15 posts, small, with Depth Marker
PL-M15	Paralight 15 posts, medium, with Depth Marker
PL-L15	Paralight 15 posts, large, with Depth Marker

PARMAX PARALIGHT REAMERS

Art no:	Product name:
	<u>2 pcs/pkg</u>
PLR-3/2	Paralight Reamers, Diam. 1,35 mm
PLR-4/2	Paralight Reamers, Diam. 1,50 mm
PLR-5/2	Paralight Reamers, Diam. 1,65 mm
PLR-6/2	Paralight Reamers, Diam. 1,80 mm

Parmax Diamonds













Super-coarse for Fast and Efficient Treatment

PARMAX CLASSIC MULTI-USE DIAMONDS

Art no:

Product name:

FG-standard, 5 pcs/packaging




201A	Parmax No. 1	
202A	Parmax No. 2	
203A	Parmax No. 3	
204A	Parmax No. 4	
205A	Parmax No. 5	
206A	Parmax No. 6	
207A	Parmax No. 7	
208A	Parmax No. 8	
209A	Parmax No. 9	
210A	Parmax No. 10	
211A	Parmax No. 11	
212A	Parmax No. 12	

PARMAX POLISHING DIAMONDS

Art no:

Product name:

FG-standard, 5 pcs/packaging

932A	Polishing diamond 378C/018	
934A	Polishing diamond C858/014	
935A	Polishing diamond C862/014	

Your Patient is Worth a Diamond!

Parmax Single-Use diamonds are useful for both the patient and the dental team.

The diamonds are manufactured according to the ISO standard.

- Disposable diamonds
- Always sharp and efficient
- No cleaning - saves time and money
- No risk of cross-contamination
- Delivered in a practical roll of blister packs.



PARMAX SINGLE-USE DIAMONDS

Art no:

Product name:

FG-standard, 50 pcs/packaging

270A	PARMAX SINGLE-USE DIAMOND 801-010M	
271A	PARMAX SINGLE-USE DIAMOND 801-016M	
274A	PARMAX SINGLE-USE DIAMOND 830-012M	
278A	PARMAX SINGLE-USE DIAMOND 835-014M	
947A	PARMAX SINGLE-USE DIAMOND 854-018M	
281A	PARMAX SINGLE-USE DIAMOND 856-016M	
282A	PARMAX SINGLE-USE DIAMOND 856-018C	
276A	PARMAX SINGLE-USE DIAMOND 858-014M	
277A	PARMAX SINGLE-USE DIAMOND 862-014M	
279A	PARMAX SINGLE-USE DIAMOND 880-014M	
273A	PARMAX SINGLE-USE DIAMOND 885-012M	
280A	S PARMAX SINGLE-USE DIAMOND 885-016M	

BL-6679-05_2025-01



A Brand of Dentatus AB, Finspångsgatan 42 163 53 Spånga
Tel: +46 (0) 8-546 509 90, e-mail: info@parmax.se
www.parmax.se